


IsoDisk®

A HEPA media filter ideal for air and gas filtration.
Suitable for use in isolators and autoclaves for the maintenance of sterile air


Features and benefits

- High efficiency particulate rating of 85% to 99.95% $\geq 0.3\mu\text{m}$
- 50mm diameter EFA 17.34cm²
- Manufactured with Borealis polypropylene
- Bi-directional air flow allowing filtration to be performed on the inlet or outlet side of the filtration device
- HEPA depth media allows for optimum flow rates
- Our tamper-printing technology allows us to incorporate your logo onto the filter housing, offering brand marketing opportunities.
- Media is ultrasonically welded.

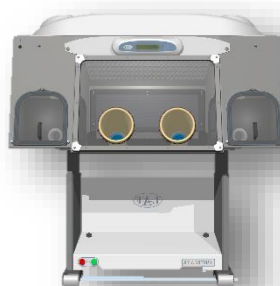

Testing available on request

- Clean pressure drop
- Water-entry-pressure
- Flow rate
- Burst integrity test
- Pressure test

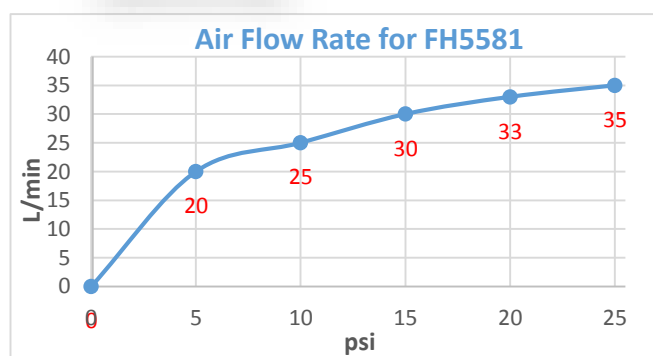

Order	Micron Size	Media	Connector	Size	Sales Unit
FH1551	0.3 μm	HEPA	1/4" stepped barb	50mm	20 filters
FH5581	0.3 μm	HEPA	3/16" straight barb	50mm	20 filters
FH25101	0.3 μm	HEPA	1/8" BSPT cut flush	50mm	20 filters
FH0371	0.3 μm	HEPA	3/16" hose barb	25mm	50 filters